


INSIGHT

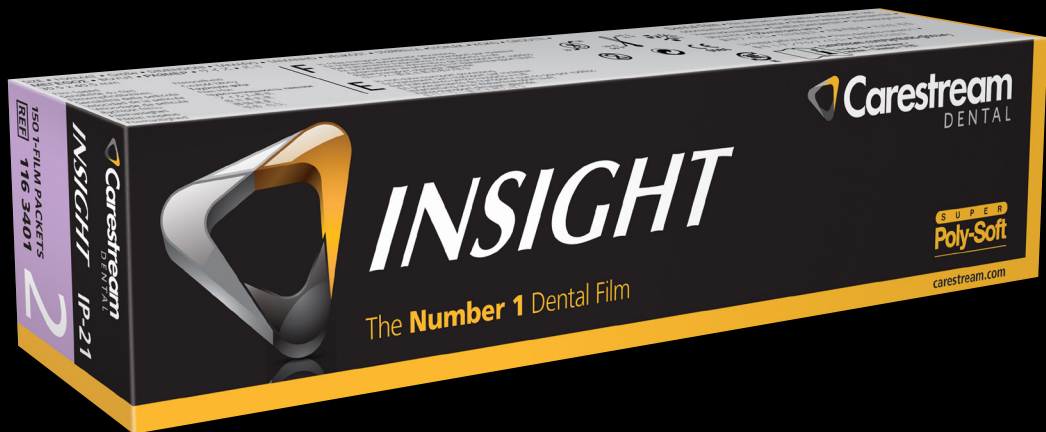

Minimize exposure. Maximize results.

The standard for intraoral film. With INSIGHT dental film, you'll get a true F-speed film that requires up to 60% less radiation per exposure.

- Delivers excellent and consistent image quality through T-Grain emulsion technology.
- Lowers patient and staff radiation exposure without sacrificing image resolution.
- Needs no continual equipment adjustments to accommodate processing conditions.
- Holds speed and contrast in variable chemistry conditions.
- Is available in all intraoral film sizes and formats, including ClinAsept barrier packets for added safety.

Exposure Guidelines

General D-speed film users		General E-speed film users	
Automatic processing	Reduce exposure time 60%*	Automatic processing	Reduce exposure time 20%*
Manual processing	Reduce exposure time 50%*	Manual processing	Reduce exposure time 20%*

*Reduce exposure time only. Do not change mA or kVp settings. To fine tune image – check exposure recommendations inside box.

For best results, use with READYMATIC, READYPRO, or GBX processing chemicals. INSIGHT dental film processes as an F-speed film when processed in roller transport automatic processors; on the high-end of E-speed range when processed manually.

Ordering Information for INSIGHT

Size	Description	Code	Films/Pkgs	REF No
Size 0	Super Poly-Soft packets	IP-01	100	867 5332
Size 0	Super Poly-Soft packets (double film)	IP-02	100	128 0619
Size 0	Super Poly-Soft packets with ClinAsept barrier	IP-01C	75	120 0328
Size 1	Paper packets	IP-11	100	112 4981
Size 1	Paper packets (double film)	IP-12	100	805 5402
Size 2	Super Poly-Soft packets	IP-21	150	116 3401
Size 2	Super Poly-Soft packets (double film)	IP-22	130	179 8628
Size 2	Super Poly-Soft packets with ClinAsept barrier	IP-21C	100	153 9931
Size 2	Super Poly-Soft packets with ClinAsept barrier (double film)	IP-22C	100	156 0390
Size 2	Paper packets	IP-21	150	176 3960
Size 2	Paper packets (double film)	IP-22	150	107 9086
Size 2	Bite-wing paper packets	IB-21	50	180 7650
Size 3	Bite-wing paper packets	IB-31	100	829 8929
Size 4	Occlusal paper packets	IO-41	25	116 9143

For more information, contact your local authorized representative for Carestream Dental products or visit: carestream.com.



Carestream Health France
1, rue Galilée
93192 NOISY-LE-GRAND CEDEX
FRANCE